

Uno Vita Red Light Therapy Cap

Instructional Manual

Red-Light-Therapy for Hair Regrowth and Decrease Pain

Thank you for purchasing our red light cap



Say hello to your red light therapy cap!

Before starting your first session, please make sure to read this manual as it may explain important operating procedures.

Notes: Each led diode inside with three chips at 1 pc 630nm, 1 pc 850nm and 1 pc 940nm. Kindly noted that near infrared led is invisible by naked eyes.

This package includes

- 1* Red light cap
- 1* Power cord
- 1* Adapter
- 1* USB Cable
- 1* User Manual


Specifications

Rating Power	75W
Actual Power	26.1W / 19.6W /14.5W
Diameter	200 MM
LED QTY	150*0.5W
Dimmable	Yes (3 levels)
Net Weight	0.18KG (for cap only)
Wavelength	1*630nm; 1*850nm; 1*940nm
EMF	0 UT
Timer	15 Minutes
Size Adjustable	Yes, with zipper
Material	Fabric
Pulse	Yes at 10 Hz
Dual Buttons	ON / OFF Dimmable Pulse Switch wavelength
Adaptor	CE ROHS UL

The red-light device uses combo red 630nm and Infrared 850nm & 940nm LEDs which is chosen to give the ideal wavelengths found in the scientific research of the last three decades.

User Instruction

1) Plug into the AC socket or powered by power bank (not included) with USB cable and the two buttons illuminate **cyan** then press the power button to turn the device on. (It comes to lowest level for brightness 50%.)

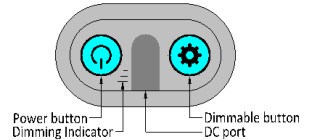

2) Distance: zero inches for best results, place your head for hair re-growth, headache and improve brain memory and brain disorder.

3) Select options:

Press the power button once quickly so the power button illuminates **green**, then:

a. Press the button once, and all LEDs (630nm 850nm 940nm) turn on.

b. Press again, and only the red LEDs (630nm) will be on, and near infrared (850nm and 940nm) off.

c. Press again, only near infrared LED (850nm and 940nm) on, and red LEDs (630nm) off. **Note that infrared light is invisible to the human eye so it may appear that the LEDs are off when in fact the infrared led is on.**

d. Press again, and all LED diodes turn off.



If you press the power button for more than 3 seconds, the power button illuminates **blue** and:



- a. Press the button once, and all LEDs (630nm 850nm 940nm) come to pulse mode.
- b. Press again, and only the red LEDs (630nm) come to pulse mode.
- c. Press again, and only the near infrared LEDs (850nm and 940nm) come to pulse mode.
- d. Press again, and all LED diodes turn off.

For **dimnable mode adjusting brightness**, please switch to the other button.



Model 1: Low level (power: 50%)

Model 2: Middle level (power: 75%)

Model 3: High level (power: 100%)

When the device is on it comes to low level already, press again, it comes to middle level, press again it comes to high level. For which levels to choose depends on you and the season. Suggest highest level for winter season and lowest level for summer season.

4) Usage Frequency: Use device 1-2 times per day for at least 4 weeks, and for convenience the device automatically shuts off after each 15-minute session. Let the device cool for a few minutes before another treatment.

5) Unplug device from wall and store in cool dry area.

6) Maintenance: For cleaning we recommend using a damp towel or wet tissue to **slightly** clean your cap. Make sure the device is disconnected from the power source, cords disconnected and turned off. Only clean the fabric portion of the cap and not the LED diodes directly. Do not submerge the cap in water. Allow the red-light cap to air dry and make sure to keep it out of direct sunlight.

Benefits backed by science :

Research shows, that the wavelengths in the spectrum of 630-950nm show optimal biotical response.

- **Decreases pain or headache**
- **Hair Re-growth**
- **Improve Brain Disorder**
- **Improve Brain Memory**

Even though this technology does not have side effects, we suggest taking the following precautions:

- ▶ **This device is water-resistant but not waterproof ! Do not submerge the wrap or its components in water.**
- ▶ **Not intended for the use by children and pregnant woman.**
- ▶ **Do not exceed 30 Min. on one treatment.**
- ▶ **Don't stare directly at the light .**
- ▶ **Contact your doctor if you :**
 - Have any medical questions.**
 - Suspect that you have had a negative reaction from using the device.**
 - Are taking drugs known to cause light sensitivity.**

Application

Here is one example of how you can use the red-light device. You can enjoy the red light at home.





Distributor: Uno Vita AS

Phone: + 47 22 09 18 80

E-mail: sales@unovita.com

Homepage: shop.unovita.com